

Mighty in Miniature ***Space-Pack & Space-Fill***


www.eighteeth.com





Magnetic Fluency Consistent Power

Milliohm-level, low resistance for rapid heating and sustained power
Magnetic design enables swift, hassle-free exchange



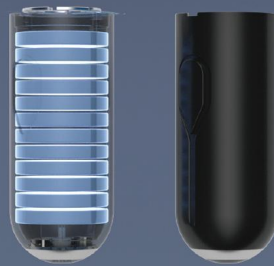


Dual Battery Configuration for All-Day Productivity

Space-Pack: Powers 800 root canal treatments per battery
Space-Fill: Powers 70 root canal treatments per battery*

Lightweight Pen-Style Design





- Slim, ergonomic body for comfortable, effortless grip
- Facilitates easy, precise handling during clinical procedures

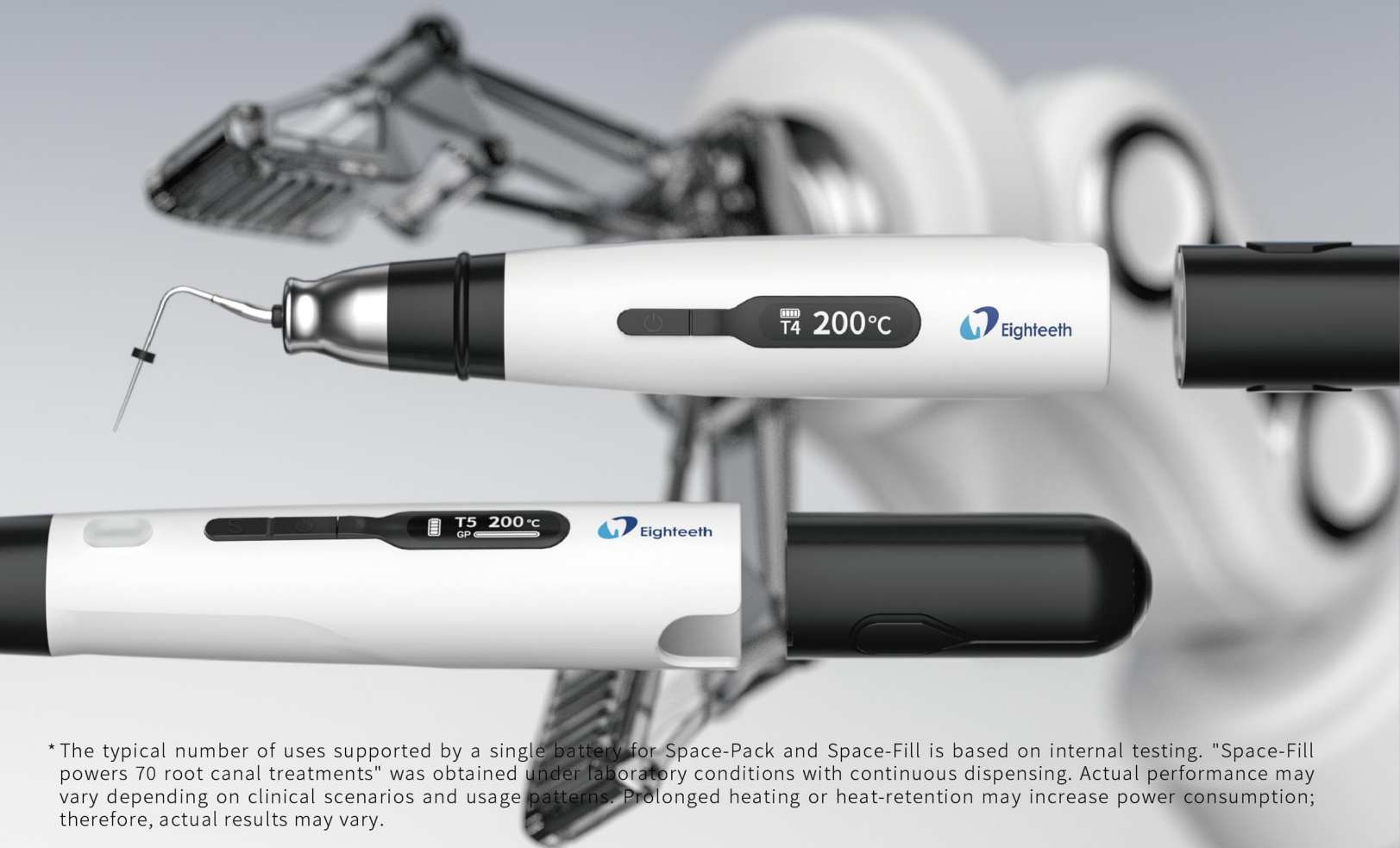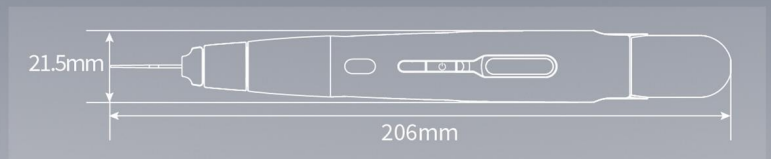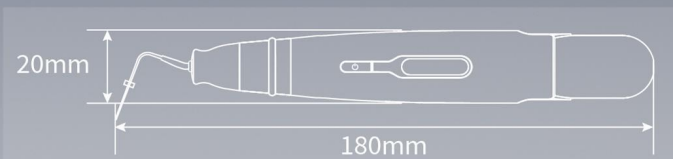

* The typical number of uses supported by a single battery for Space-Pack and Space-Fill is based on internal testing. "Space-Fill powers 70 root canal treatments" was obtained under laboratory conditions with continuous dispensing. Actual performance may vary depending on clinical scenarios and usage patterns. Prolonged heating or heat-retention may increase power consumption; therefore, actual results may vary.

Comprehensive Working Tips for Precision Dentistry


Heat Tip

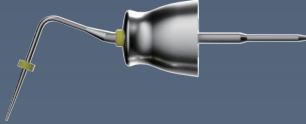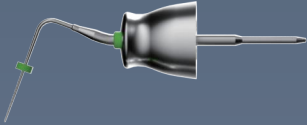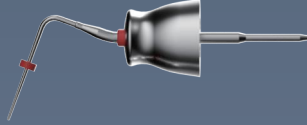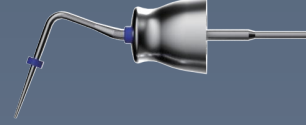Standard *(one for each)	Tip 40 ϕ 0.4mm	Tip 45 ϕ 0.45mm	Tip 50 ϕ 0.5mm
			
Optional	Tip 35 ϕ 0.35mm	Tip 55 ϕ 0.55mm	Tip 60 ϕ 0.6mm
			

Fill Tip

Standard *(one for each)	23Ga	25Ga	
			

* The specification data of Space-Pack and Space-Fill is based on internal testing. Due to differences in actual usage scenarios, operating methods, and other factors, the specific values may vary. Please refer to actual usage conditions for accuracy.

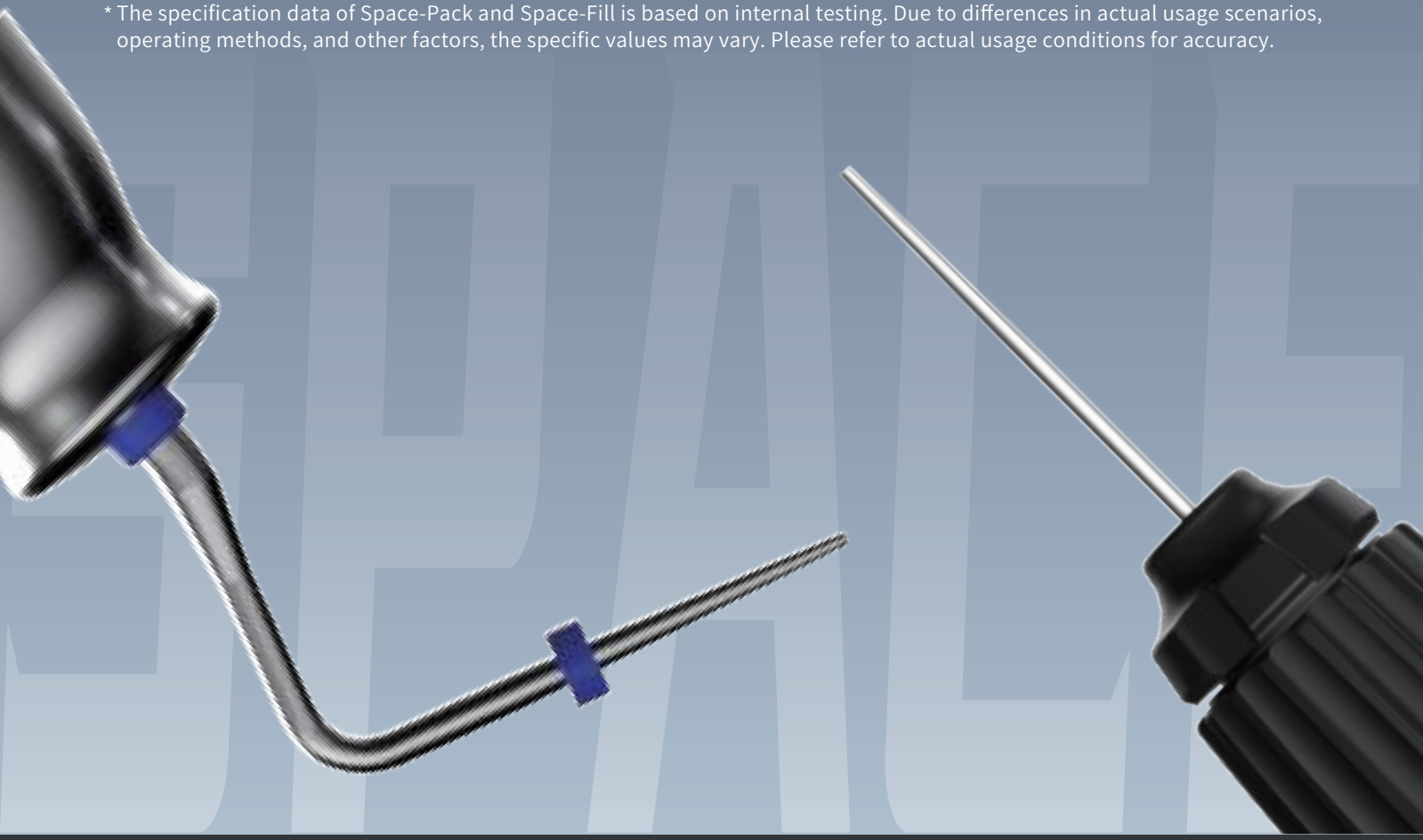

Space-Fill

360° Rotatable Autoclavable Silver Needle



- Innovative coupling technology prevents leakage and overflow
- 999 pure silver maintain thermoplasticity for smooth flow
- Pre-Bent, extra-long design for convenience and visibility

Electric Motorized GP Flow



- Instant start/stop, retraction prevents leakage
- Smooth, bubble-free extrusion for dense & precise filling

Five Memory Settings (T1-T5)



- Quick access to personalized configurations

Precision Temperature Control



- 3 adjustable filling speeds: low, medium, high
- Temperature range: 100°C - 200°C

LED Display Screen



- Real-Time monitoring of gutta-percha remaining






Instantaneous Heating

Heating to 180°C takes only
10 seconds



Leakage-Proof Outlet Design

Innovative extrusion outlet
structure effectively prevents
leakage and simplifies the
cleaning process

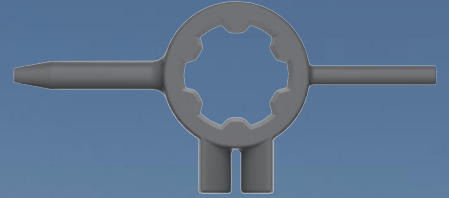

Standard Multi-Function Wrench

- Outlet Residue Removal & Inner
Tip Cleaning
- Needle Angle Rotation
- Bending adjustments without breakage

Innovative Two Filling Modes for Diverse Clinical Needs

- Standard Filling Mode
- Inject Control Mode with five selectable settings

Standard Insulating Sleeve for Enhanced Safety **Comprehensive Fill Tip Selection**



* "Innovative Two Filling Modes" refers to the addition of a inject control mode to the standard filling mode of the Space-Fill.

Space-Pack



Upgraded Heating Tip

- Targeted 4-5mm tip end heating for precision cutting
- High toughness guarantees enduring reliability



360° Ring Button

- Sensitive control, precise left-right handed operation
- Complete with audible feedback for each press



Five Memory Settings (T1-T5) for Customized Operation



Real-time Temperature Display



Automatic Power-Off Protection

- Adjustable heating time: 3-10 seconds
- Customizable power-off delay for specific procedures



Smart Temperature Control System

- Adaptive compensation algorithm ensures consistent & stable temperatures, and enhances cutting efficiency in challenging conditions including moist environments and complex loads
- Wide adjustable temperature range: 90°C to 250°C for diverse filling materials

Heating Faster, Waiting Less

- Heats up to 200°C in just 0.07 seconds

Diverse Heat Tip Options

- Specialized tips cater to versatile clinical needs
- Tailored heating profiles optimize performance for each application



* The heating time data is based on internal testing using the Space-Pack working tip with a diameter of $\Phi 0.35\text{mm}$. Due to differences in actual usage scenarios, operating methods, and other factors, the specific values may vary. Please refer to actual usage conditions for accuracy.

Specification

Space-Pack

Heating Time	0.07 seconds for 200 °C (Working Tip: 35)
Charging Time	1 hour to complete charge
Operating Time	Up to 800 treatments per full charge
Weight (Include Tip)	84g
Length (Include Tip)	180mm

Space-Fill

Heating Time	10 seconds for 180°C(23Ga/25Ga)
Charging Time	1 hour to complete charge
Operating Time	Up to 70 treatments per full charge
Weight (Include Tip)	129g
Length (Include Tip)	206mm

